



**BELDICO** ... *We simply love babies*



## The disposable feeding bottles

*'The standard in medical disposable feeding bottles'*



## INTENTIONS OF USE

- We have a full range of feeding bottles for children, which fulfil the needs in **neonate, maternity and paediatric departments.**

- All our feeding bottles are ready to use.

- The complete range of disposable feeding bottles and teats is available **in various configurations** which are **adapted to each specific hospital need.**

- The Beldico disposable feeding bottles are **used for the collection and storage of mother milk**, as well as for the **preparation of formula milk.** **The milk can also be pasteurised** in our feeding bottles.

- Our full range of disposable teats, **made of injected (latex-free) polymer**, is adapted for **children of all age** and **fits all the Beldico feeding bottles.**



**SAFE**

The Beldico feeding bottles are produced in food-grade quality polypropylene. The bottles are in accordance with the European Directives for plastic containers for foodstuff.

The Beldico feeding bottles are leak proof. The calibrated graduation offers a reliable accuracy.

The Beldico feeding bottles are manufactured, assembled and packed under clean and controlled production conditions. (ISO 13485 and 9001)

The thermoplastic teats are produced in a neutral material that does not contain toxic derivatives.



**BELDICO**

... **WE**  
**SIMPLY**  
**LOVE**  
**BABIES**





## PRATICAL USE

The Beldico disposable feeding bottles are available in 50, 120 and 240 ml.

Our feeding bottles can be frozen for storage to -30°C and can be microwave heated.

Our feeding bottles do also fit in the standard contact heaters and transport baskets.

As all the teats fit all the bottles sizes, there is always a complete bottle that meets the child's requirements.




## HYGIENIC

As the feeding bottles, caps and teats are ready to use, there is no risk of contamination due to cleaning and decontamination treatments in the hospital.

The Beldico feeding bottles are supplied in peel-open pouches and are gamma treated or EO sterilized.

The teat protectors assure an aseptic cover.

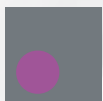

## ECOLOGICAL

The plastic of the bottles, rings, caps and teats can be recycled if wanted.

The bottles, rings, caps and teats can be incinerated without any toxic emanations.



## ECONOMICAL




As the feeding bottles, teats and caps are available as single use, disposable products, there is no need to have cleaning and sterilisation procedures, room, nor equipment.

The products are immediately available for distribution of the required food.



## OUR TAMPER-PROOF CAPS

Our tamper-proof caps guarantee a closed system and perfect sealing during transport, storage and pasteurization.

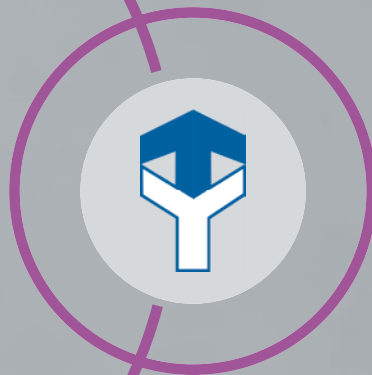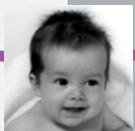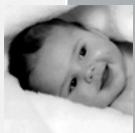

## OUR ENTERAL FEEDING CAP

With our specially designed enteral feeding cap, the Beldico feeding bottles can also function as an enteral nutrition bottle.

A special hanger is available to fix the feeding bottle on an infusion standard.

**NEW!** A 500 ml bottle with double graduation is available for enteral feeding preparations. The bottle has reversed graduation and an eyelet for easy hanging.



## OUR TEATS

The Beldico (latex-free) teats are produced in a synthetic thermoplastic compound, commonly used in the medical and food industry.

For premature babies, we have developed a special 3-speed teat, which guarantees an easy feeding of all nutrients (pink ring).

For newborn infants, we have a 1-speed teat with an anatomical shape (red ring).

For newborns infants and babies in paediatric, we have created a 3-speed teat, which guarantees a correct flow rate following the needs of the child (blue ring).

For newborns and small babies in paediatric, who temporarily cannot receive breast-feeding, but who will afterwards go or return to breast-feeding, we have developed a specially shaped anatomical teat (white ring). This teat prevents nipple confusion.

## The Beldico teats...



*A good and comfortable hold in the baby's mouth.*



*Completely odour-free and flavourless.*



## OUR TEAT PROTECTORS

The Beldico teat protectors have a unique and reliable click system, which ensures a perfect seal for the teat/ring. Transport and handling at the hospital can take place without any risk. Moreover, aseptic handling is fully assured.

The graduated teat protectors also function as a milk cup to be used in the neonatology and maternity ward.



**BELDICO S.A.**

rue Andre Feher 1  
Zoning de Aye  
6900 Marche-en-Famenne  
BELGIUM

Tel.: +32 (0)84 31 37 80

Fax: +31 (0)84 31 37 87

E-mail: [info@beldico.be](mailto:info@beldico.be)

Website: [www.beldico.be](http://www.beldico.be)

*Distributed by:*

